



Dose-volume histogram — Lung_Patient_3

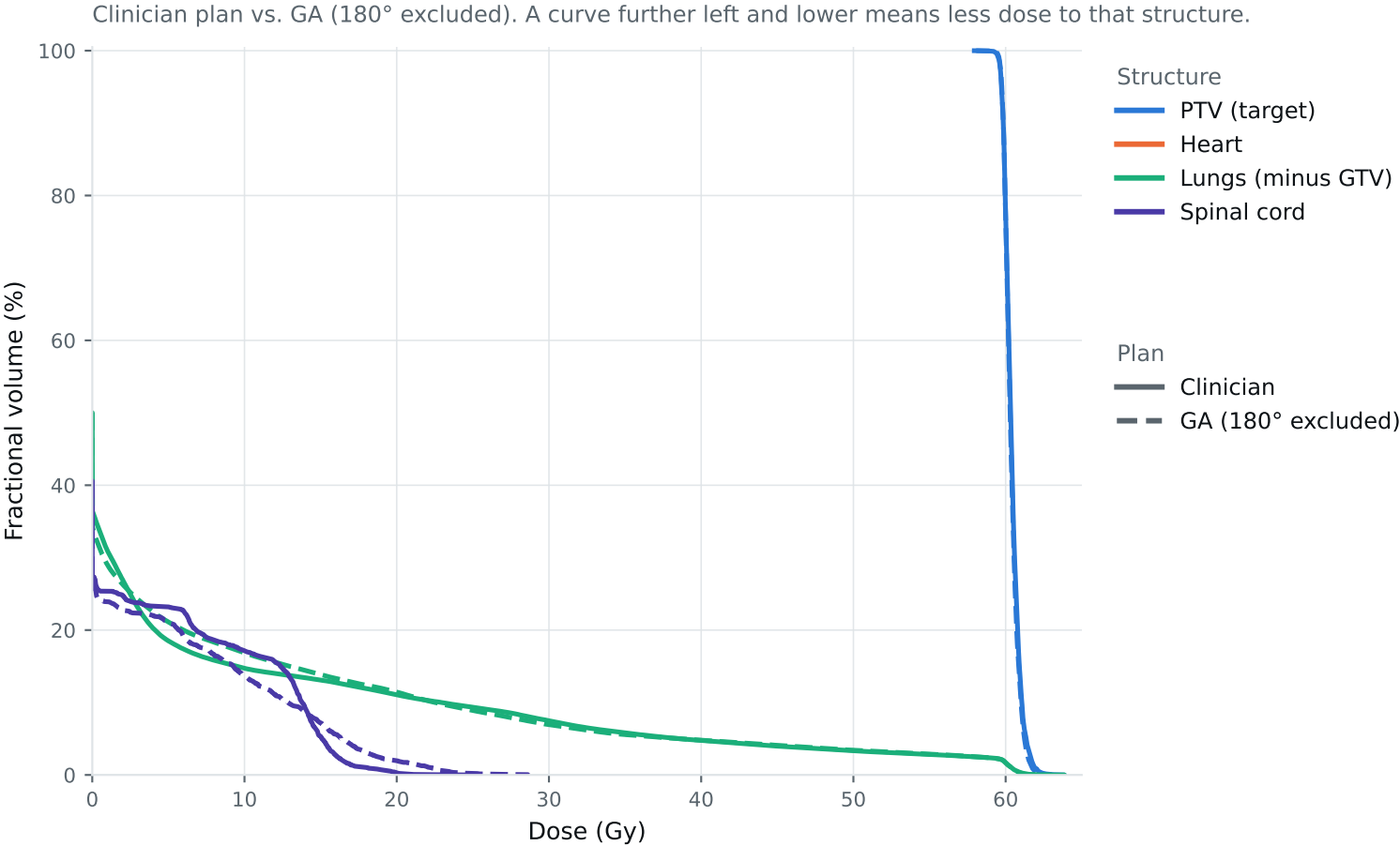

Curves further left are better for organs; the PTV curve should stay high, then drop sharply. Target curves should overlap; organ curves show the tradeoffs. Esophagus is on page 2.

DVH by structure — Lung_Patient_3

Same curves, one structure per panel, all plans, each annotated with that structure's protocol metrics.
Bold = best plan; grey = all plans are within 0.01 of each other, so the row carries no signal. Amber misses a goal, red exceeds a limit.

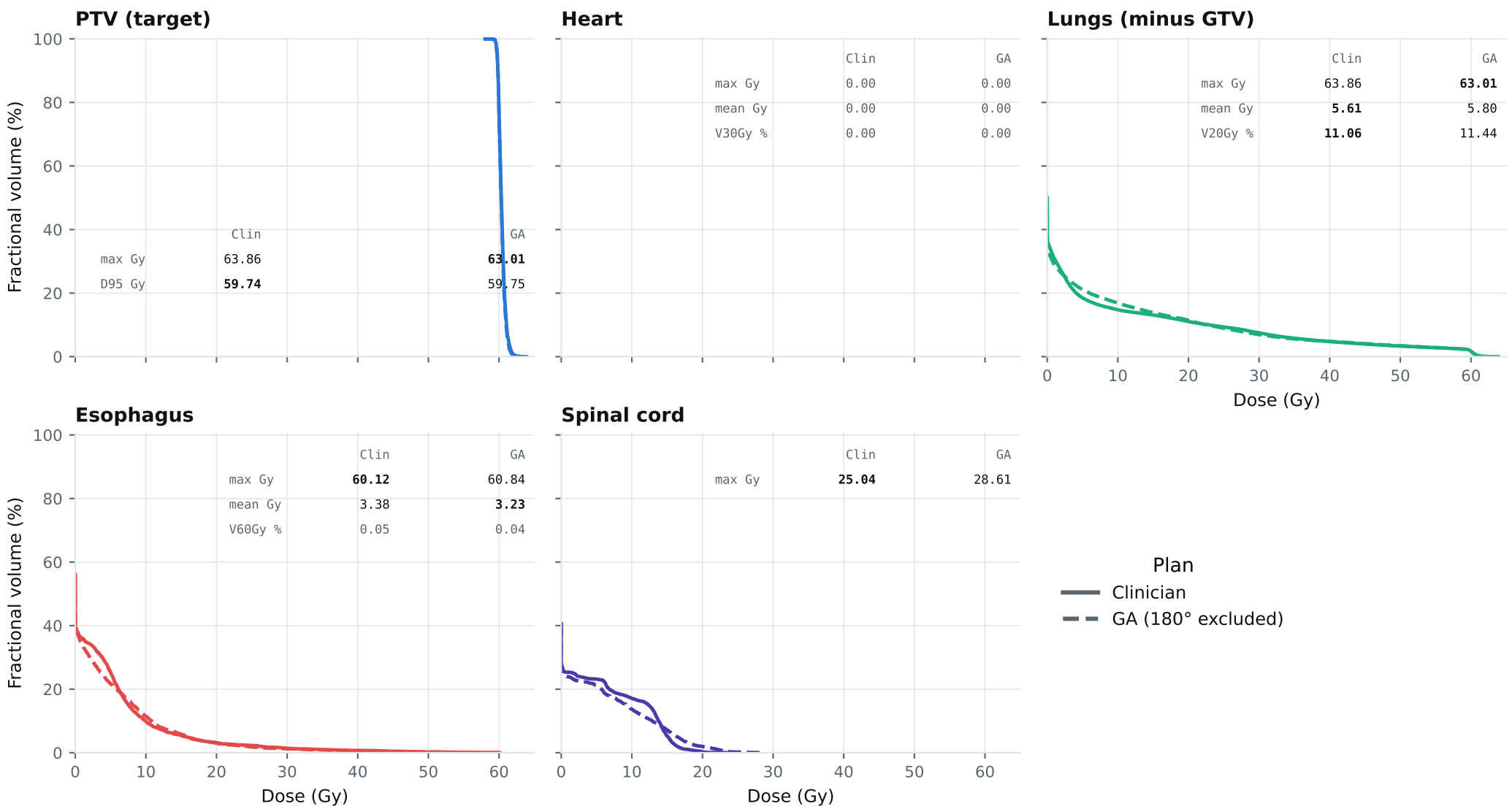

Clinical criteria — Lung_Patient_3

Lung_2Gy_30Fx protocol. Every criterion is a ceiling, so lower is better; bold marks the best plan.

Criterion	Limit	Goal	Clinician	GA 180° excluded	
GTV max	69	66	61.89	61.87	Gy
PTV max	69	66	63.86	63.01	Gy
ESOPHAGUS max	66	—	60.12	60.84	Gy
ESOPHAGUS mean	34	21	3.38	3.23	Gy
ESOPHAGUS V60Gy	17	—	0.05	0.04	%
HEART max	66	—	0.00	0.00	Gy
HEART mean	27	20	0.00	0.00	Gy
HEART V30Gy	50	48	0.00	0.00	%
LUNG_L max	66	—	63.86	63.01	Gy
LUNG_R max	66	—	28.27	33.17	Gy
CORD max	50	48	25.04	28.61	Gy
SKIN max	60	—	45.35	44.95	Gy
LUNGS_NOT_GTV max	66	—	63.86	63.01	Gy
LUNGS_NOT_GTV mean	21	20	5.61	5.80	Gy
LUNGS_NOT_GTV V20Gy	37	—	11.06	11.44	%
PTV D95 (target coverage)			59.74	59.75	
Objective value (search signal)			53.22	52.50	
Solve time			25 s	23 s	
Beam angles (gantry°)	Clinician	0, 30, 60, 90, 120, 150, 175			
	GA	0, 45, 75, 150, 165, 210, 300			

Grey rows: all plans agree within 0.01 — heart, lung and LUNGS_NOT_GTV max sit at 66 Gy because those structures abut the target. Amber misses the goal; red exceeds the limit.

Source: Lung_Patient_3_clinical_metrics_full_resolution.json, full resolution. The curves and table come from the same solve.