



Dose-volume histogram — Lung_Patient_5

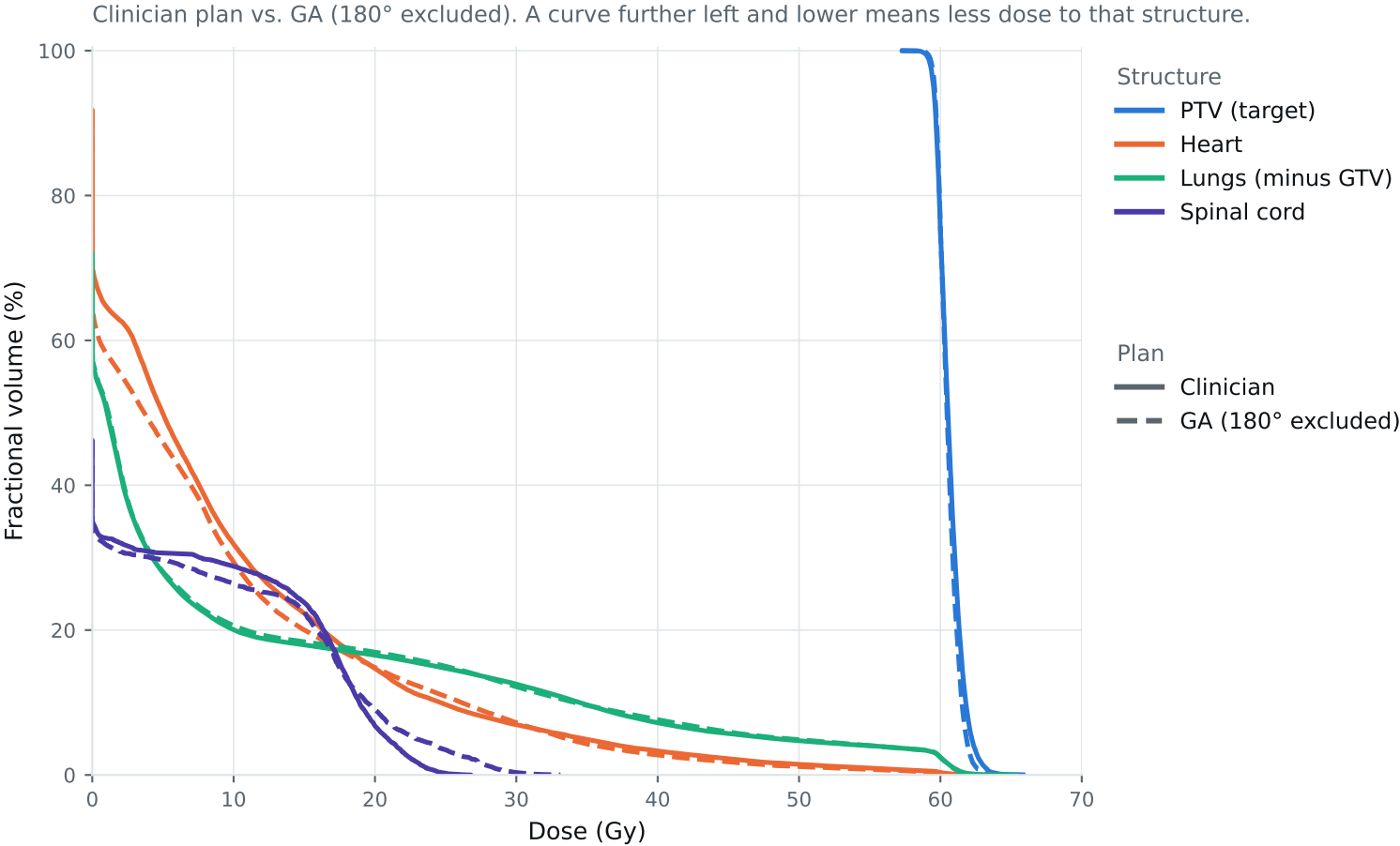

Curves further left are better for organs; the PTV curve should stay high, then drop sharply. Target curves should overlap; organ curves show the tradeoffs. Esophagus is on page 2.

DVH by structure — Lung_Patient_5

Same curves, one structure per panel, all plans, each annotated with that structure's protocol metrics.
Bold = best plan; grey = all plans are within 0.01 of each other, so the row carries no signal. Amber misses a goal, red exceeds a limit.

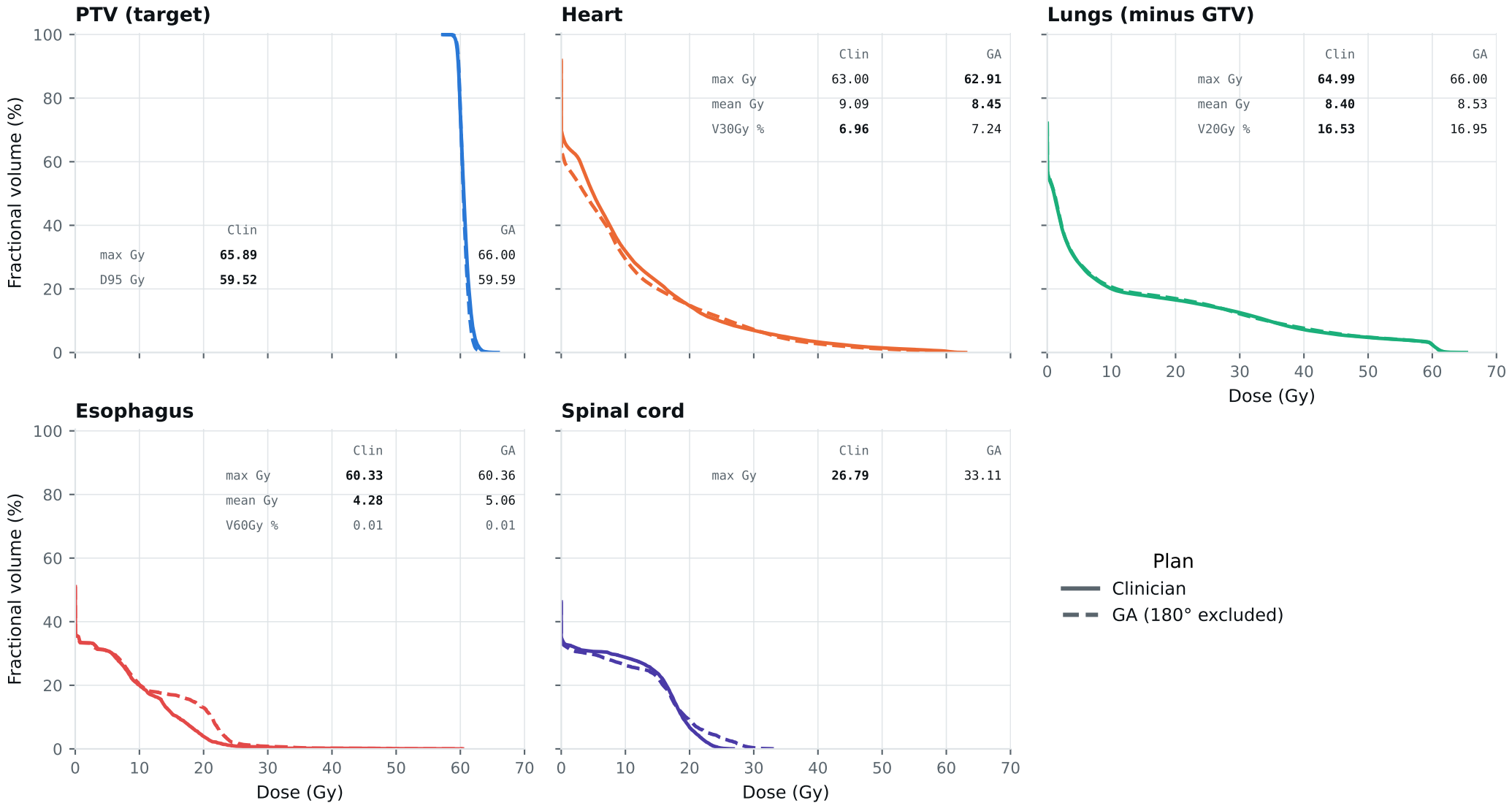

Clinical criteria — Lung_Patient_5

Lung_2Gy_30Fx protocol. Every criterion is a ceiling, so lower is better; bold marks the best plan.

Criterion	Limit	Goal	Clinician	GA 180° excluded	
GTV max	69	66	64.06	63.82	Gy
PTV max	69	66	65.89	66.00	Gy
ESOPHAGUS max	66	—	60.33	60.36	Gy
ESOPHAGUS mean	34	21	4.28	5.06	Gy
ESOPHAGUS V60Gy	17	—	0.01	0.01	%
HEART max	66	—	63.00	62.91	Gy
HEART mean	27	20	9.09	8.45	Gy
HEART V30Gy	50	48	6.96	7.24	%
LUNG_L max	66	—	64.99	66.00	Gy
LUNG_R max	66	—	40.00	38.73	Gy
CORD max	50	48	26.79	33.11	Gy
SKIN max	60	—	45.00	43.72	Gy
LUNGS_NOT_GTV max	66	—	64.99	66.00	Gy
LUNGS_NOT_GTV mean	21	20	8.40	8.53	Gy
LUNGS_NOT_GTV V20Gy	37	—	16.53	16.95	%
PTV D95 (target coverage)			59.52	59.59	
Objective value (search signal)			79.98	74.78	
Solve time			48 s	48 s	
Beam angles (gantry°)	Clinician	0, 30, 60, 90, 120, 150, 175			
	GA	15, 60, 90, 120, 175, 210, 345			

Grey rows: all plans agree within 0.01 — heart, lung and LUNGS_NOT_GTV max sit at 66 Gy because those structures abut the target. Amber misses the goal; red exceeds the limit.

Source: Lung_Patient_5_clinical_metrics_full_resolution.json, full resolution. The curves and table come from the same solve.