



Dose-volume histogram — Lung_Patient_6

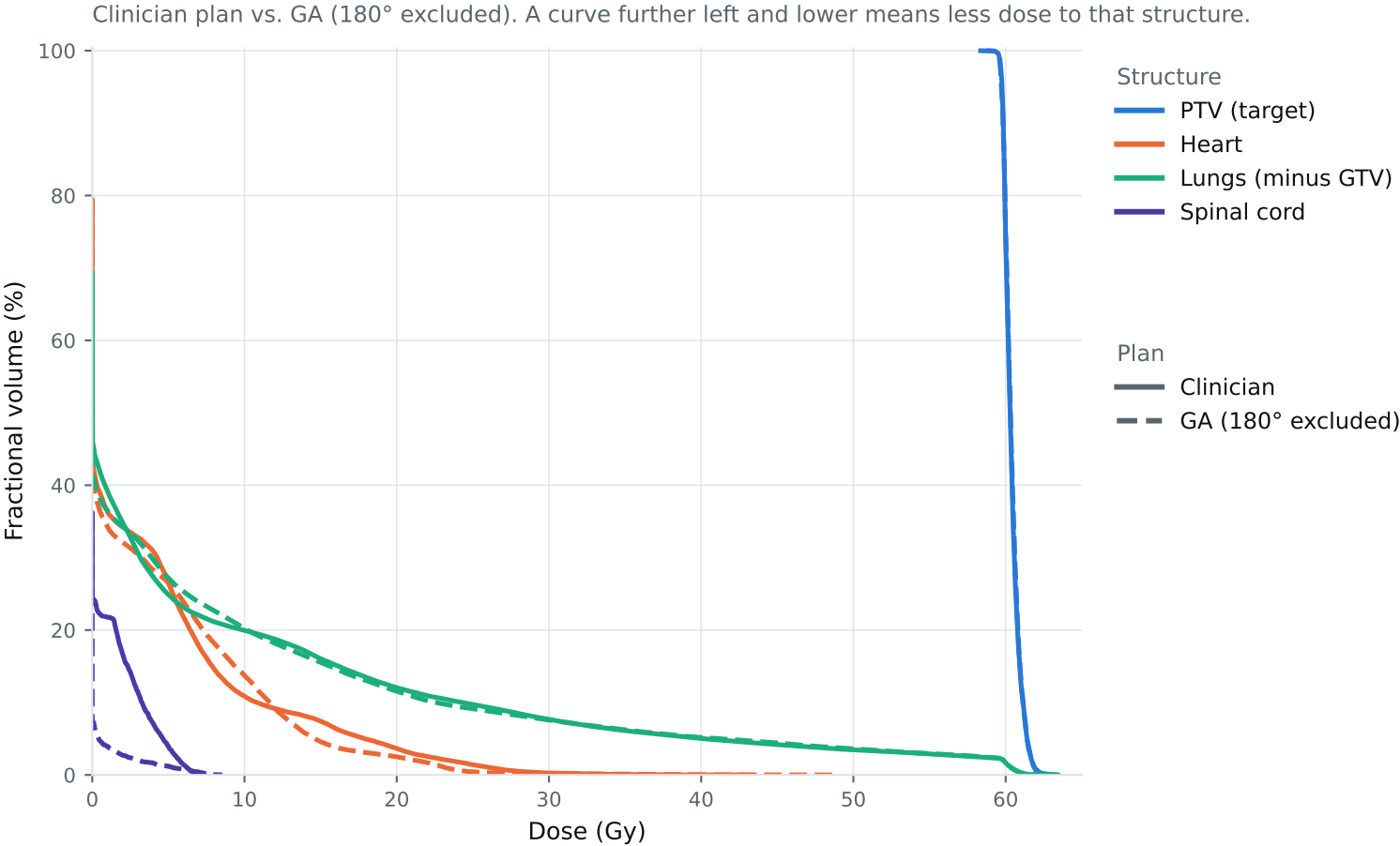

Curves further left are better for organs; the PTV curve should stay high, then drop sharply. Target curves should overlap; organ curves show the tradeoffs. Esophagus is on page 2.

DVH by structure — Lung_Patient_6

Same curves, one structure per panel, all plans, each annotated with that structure's protocol metrics.
Bold = best plan; grey = all plans are within 0.01 of each other, so the row carries no signal. Amber misses a goal, red exceeds a limit.

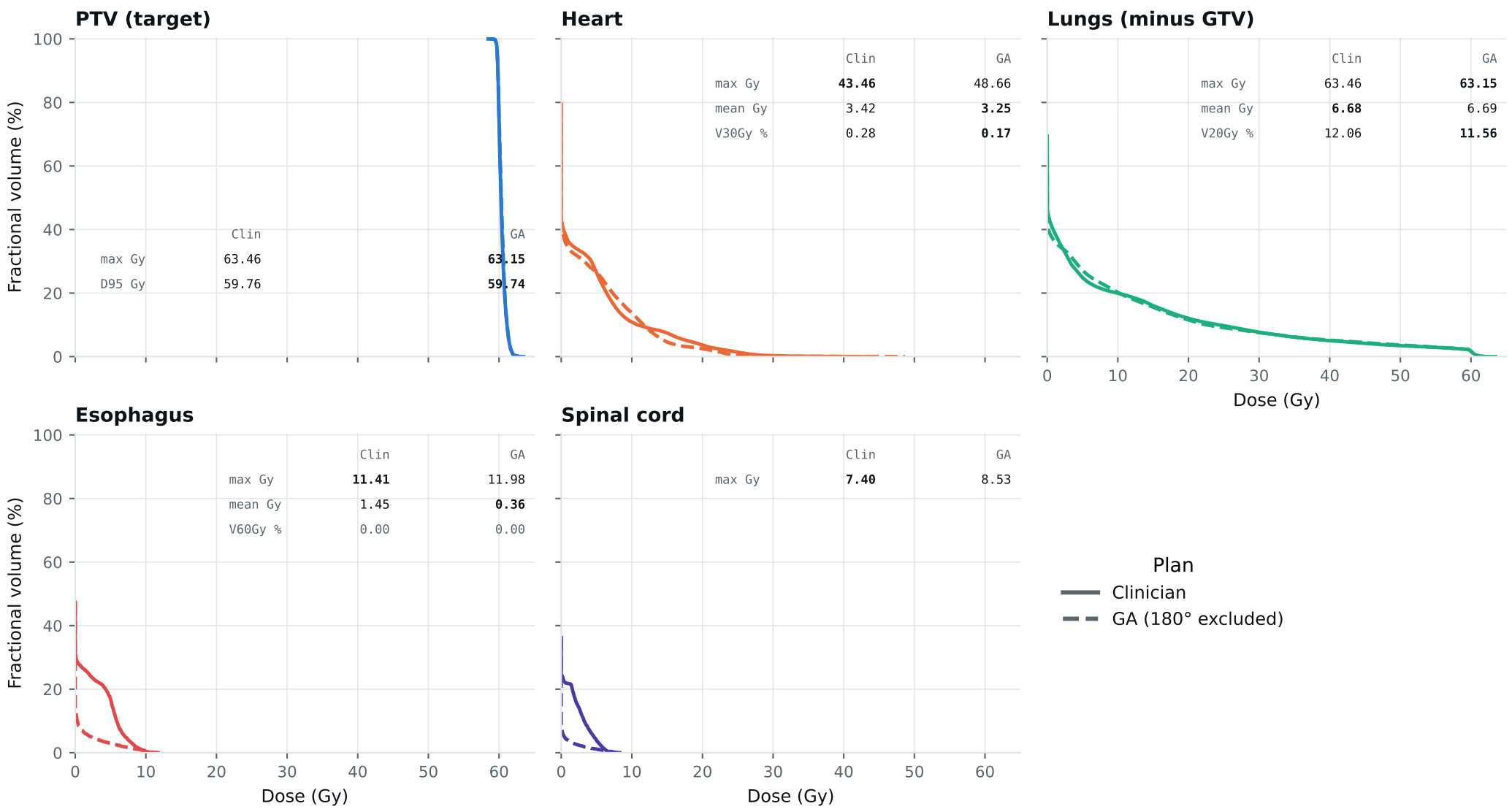

Clinical criteria — Lung_Patient_6

Lung_2Gy_30Fx protocol. Every criterion is a ceiling, so lower is better; bold marks the best plan.

Criterion	Limit	Goal	Clinician	GA 180° excluded	
GTV max	69	66	61.89	61.86	Gy
PTV max	69	66	63.46	63.15	Gy
ESOPHAGUS max	66	—	11.41	11.98	Gy
ESOPHAGUS mean	34	21	1.45	0.36	Gy
ESOPHAGUS V60Gy	17	—	0.00	0.00	%
HEART max	66	—	43.46	48.66	Gy
HEART mean	27	20	3.42	3.25	Gy
HEART V30Gy	50	48	0.28	0.17	%
LUNG_L max	66	—	11.45	26.36	Gy
LUNG_R max	66	—	63.46	63.15	Gy
CORD max	50	48	7.40	8.53	Gy
SKIN max	60	—	48.11	45.65	Gy
LUNGS_NOT_GTV max	66	—	63.46	63.15	Gy
LUNGS_NOT_GTV mean	21	20	6.68	6.69	Gy
LUNGS_NOT_GTV V20Gy	37	—	12.06	11.56	%
PTV D95 (target coverage)			59.76	59.74	
Objective value (search signal)			44.28	44.35	
Solve time			21 s	20 s	
Beam angles (gantry°)	Clinician	0, 185, 210, 240, 270, 300, 330			
	GA	15, 30, 90, 185, 240, 285, 345			

Grey rows: all plans agree within 0.01 — heart, lung and LUNGS_NOT_GTV max sit at 66 Gy because those structures abut the target. Amber misses the goal; red exceeds the limit.

Source: Lung_Patient_6_clinical_metrics_full_resolution.json, full resolution. The curves and table come from the same solve.